

DENTAPEN®




DENTAL ANESTHESIA
Painless Effortless Cordless

FOR YOUR PATIENTS


“The number 1 criteria in choosing a dentist is:
A dentist who can provide painless anesthesia.”

How dentists are judged by patients, St Georges J, Dent Today, 2004.

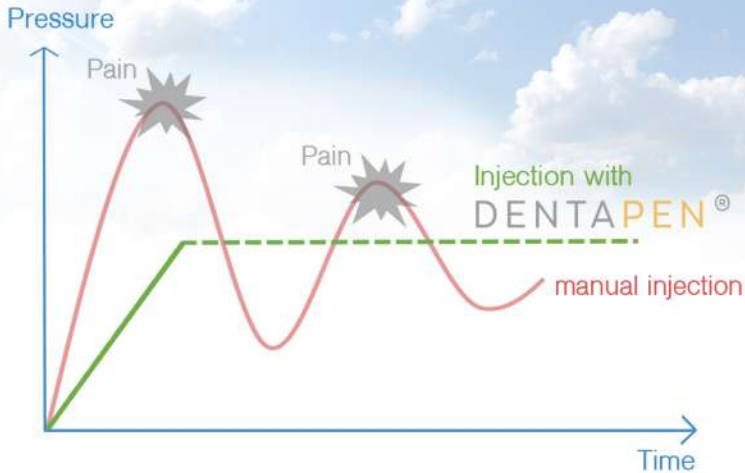

“The results of this review indicate that using C-CLAD
resulted in less pain and more effective anesthesia.”

*Computer-controlled local anesthetic delivery for painless anesthesia: a literature review,
Eun-Jung Kwak, Nan-sim Pang, Jin-Hyung Cjo, Bock Young Jung, Kee-Deog Kim and Wonse Park
Dental Anesth Pain Med. 2016 Jun; 16(2):81-88.*

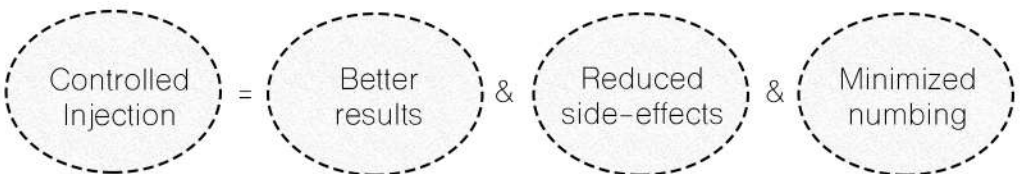

DENTAPEN®

ROTATING BARREL
for easy bevel placement



BATTERY
up to 100 patients

2-YEARS WARRANTY


AUTOCLAVABLE
low operating costs

Complete Dentapen® Kit:

	Ref.
Metric hub	DK01.2/M
Imperial hub	DK01.2/ I

Optional:

Protective sleeve x 250	PS01
-------------------------------	------

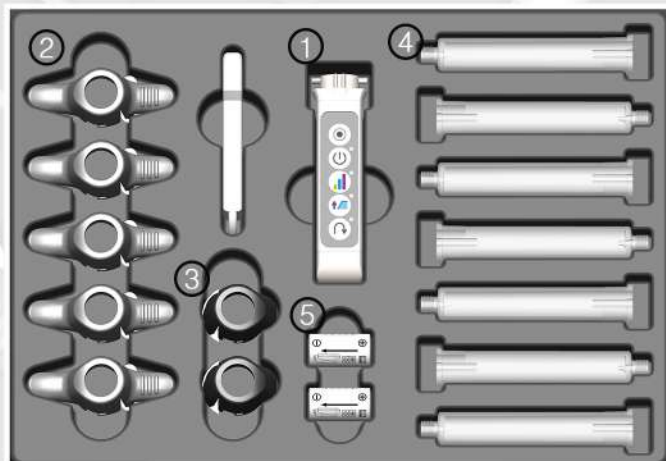

- 1 Dentapen® x1
- 2 Finger grip x5
- 3 Finger grip for nerve block x2
- 4 Cartridge holders 1.7ml metric or imperial hub x7
- 5 Battery x3

Includes:
25 protective sleeves

Designed for EFFICIENCY



Autoclavable accessories

Universal: all techniques, metric & imperial needles,

Two-way activation : index or thumb activation for more comfort

Battery powered

Compact and light (50g. with battery)

Rotating barrel for easy bevel placement

Easy Set-up



Power On/Off
Battery charge level

☀ full capacity 🌟 20% 🔴 10%



☀ Slow 90 s/mL
🌟 Medium 60 s/mL
🔴 Fast 30 s/mL



🔴 Intraligamentary mode
🌟 Ramp-up mode = anesthesia flow increases gradually for totally painless procedure



Active aspiration
Rewind back to initial position

FOR YOU

Cutting-edge C-CLAD* technology for advanced dental anesthesia

- ✓ Controlled initial force and constant flow for smoother injections:
 - Increased anesthesia success rate
 - Saves time and stress
- ✓ Computerized pressure feedback sensor:
 - Adjusts the injection speed depending on pressure of the tissue.
- ✓ Other benefits:
 - Allows challenging injections
 - Performs aspirations for a safe procedure
 - Decreases risk of work related injuries: tendinitis, arthritis...



LESS

MORE



PAIN



FEAR



NUMBING



CHILDREN
FRIENDLY



SATISFIED
PATIENTS



REFERRAL

*computer-controlled local anesthesia device



www.dentapen.ch



@dentapen_injector



Juvaplus



Juvaplus

I N J E C T I O N S Y S T E M S

J U V A +

Designed and manufactured in Neuchatel - Switzerland
Sales: Av. de Miremont, 12 - 1206 Geneve - Switzerland
info@dentapen.ch - Tel. +41 22 588 66 89