

# DENTAPEN®




**DENTAL ANESTHESIA**  
Painless Effortless Cordless

# FOR YOUR PATIENTS

“The number 1 criteria in choosing a dentist is:  
A dentist who can provide painless anesthesia.”

*How dentists are judged by patients, St Georges J, Dent Today, 2004.*

Pressure  
(N/mm<sup>2</sup>)

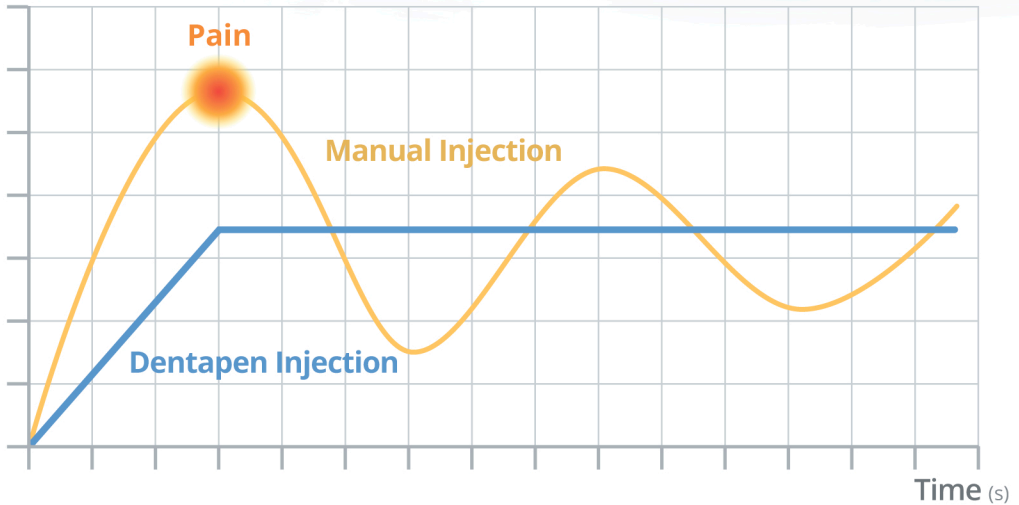

Controlled  
Injection

=

Increased  
Success Rate

+

Reduced  
side-effects

+

Minimized  
numbing

# DENTAPEN®

ROTATING BARREL  
for easy bevel placement



AUTOCLAVABLE  
low operating costs



BATTERY  
up to 100 patients


2-YEARS WARRANTY

Complete Dentapen® Kit:

Metric hub .....	Ref. DK01.2/M
Imperial hub .....	DK01.2/I

Optional highly recommended:

Protective sleeve x 250 .....	PS01
-------------------------------	------

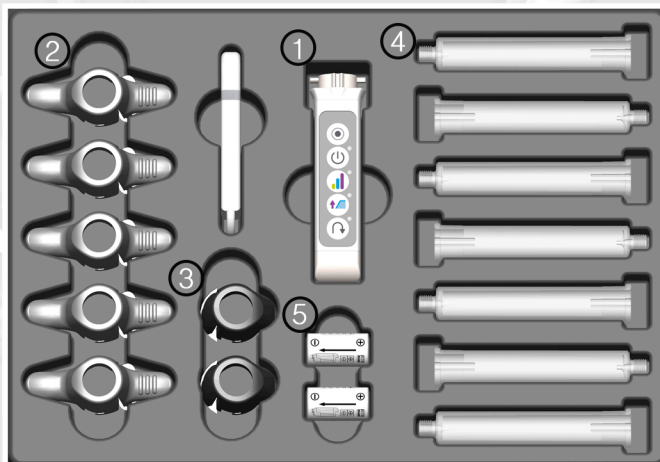

- 1 Dentapen® x1
- 2 Finger grip x5
- 3 Finger grip for nerve block x2
- 4 Cartridge holders 1.7ml metric or imperial hub x7
- 5 Battery x3

Includes:  
25 protective sleeves

# Designed for EFFICIENCY

- ✓ Autoclavable accessories
- ✓ Universal: all techniques
- ✓ Two-way activation and grip: syringe or pen like
- ✓ Battery powered
- ✓ Compact and light: only 50 grams
- ✓ Active aspiration



## Easy Set-up



- ☀ Slow 90 s/mL
- ☀ Medium 60 s/mL
- ☀ Fast 30 s/mL



- ☀ Intraligamentary mode
- ☀ Ramp-up mode = anesthesia flow increases gradually for totally painless procedure



Rewind back to initial position



# FOR YOU

## Cutting-edge C-CLAD\* technology for advanced dental anesthesia



Controlled initial force and constant flow for smoother injections:  
- Increases anesthesia success rate and saves time



Computerized pressure feedback sensor:  
- Adjusts the injection speed depending on pressure of the tissue



Other benefits:  
- Allows challenging injections: intraligamentary  
- Performs aspirations for a safe procedure  
- Decreases risk of work related injuries: tendinitis, arthritis  
- Manages patient anxiety



## LESS



PAIN



ANXIETY



NUMBING

## MORE



CHILDREN  
FRIENDLY



SATISFIED  
PATIENTS



REFERRAL

\*computer-controlled local anesthesia device



[www.dentapen.ch](http://www.dentapen.ch)



@dentapen\_injector



Juvaplust



Dentapen

R O B O T I C   S Y R I N G E S

J U V A +

Designed and manufactured in Neuchâtel - Switzerland

Sales: Av. de Miremont 12, Geneva - Switzerland

[info@dentapen.ch](mailto:info@dentapen.ch) - Tel. +41 22 588 66 89